

Oxygen Concentrator MQ50 VC5 O2

The very quiet and innovative Oxygen Concentrator

Benefits of the MQ50 Oxygen Concentrator



Exemption
from leakage
current test



Very quiet



Very
light-weight



Oxygen
24/7



Functional
display

Actual
volume level:
< 36 dB(A) *

Actual energy
consumption:
< 286 Wh **



German HMV-Nr.:
14.24.06.0028

Alarm functions

The MQ50 has the following visual and/or audible alarm functions:

- Low oxygen concentration
- Pressure deviation
- Power failure
- Low flow
- Device malfunction

Technical Data

Flow rate:	1-5 l/min
Adjustment:	Continuously variable with integrated flowmeter
Flow mode:	Constant
O ₂ Concentration:	87% - 96%
Output pressure:	0,05 +/- 0,005 MPa
Power consumption:	< 320 VA
Noise level:	< 36 dB(A)
Operating temperature:	10-37 °C
Weight:	15,5 Kg
Dimensions:	39 x 23 x 60 cm
Transport aid:	4 casters and carrying handle

Scope of delivery

- 1 x Oxygen Concentrator MQ50 VC5 O2
- 1 x Humidifier
- 2 x Nasal Cannulas
- 1 x Sponge Mesh Filter
- 1 x User manual
- 1 x PVC Tube for humidifier
- 1 x Fire Damper



Detailed front view

*Sound pressure level – A

**At setting 1 - 5 l/min



VIBRACARE GmbH
INNOVATIVE HEALTHCARE
SOLUTIONS

View of Inlet Filter



Filter access flap for Inlet Filter




Inlet Filter

View of Sponge Mesh Filter



Window/flap for Sponge Mesh Filter

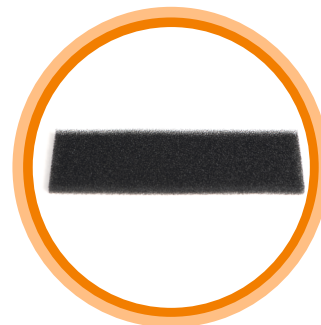


Sponge Mesh Filter

View inside



Easy access to all components
by simply opening the housing.

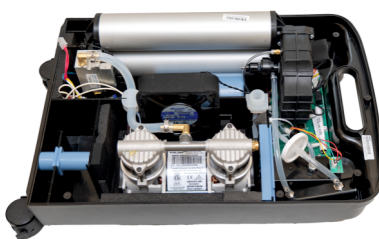

**Easy to service
and maintain**